

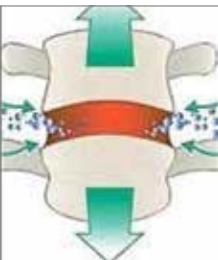
THERAMASTER

COMPONENTS



- ① Control Panel
- ② Cervical decompression
- ③ Cervical angle adjustment
- ④ Cervical one side
- ⑤ Height control
- ⑥ Emergency switch for patient
- ⑦ Up & down mobilization
- ⑧ Sacral mobilization
- ⑨ Pelvic mobilization
- ⑩ Lumbar angle adjustment
- ⑪ Lumbar decompression
- ⑫ Lumbar one side
- ⑬ Infrared / Up
- ⑭ Infrared / Down

THERAMASTER



CORE 1
• DISC DECOMPRESSION



CORE 2
• SPINE CORRECTION



CORE 3
• SPINE MOBILIZATION



CORE 4
• INTELLIGENT PROGRAM



Manufactured by Shinhwa medical
THERAMASTER
MASTER OF SPINAL THERAPY



SHINHWAMEDICAL

SHINHWAMEDICAL

SHINHWAMEDICAL

SHINHWAMEDICAL

SHINHWAMEDICAL

SHINHWAMEDICAL

SHINHWAMEDICAL

SHINHWAMEDICAL

SHINHWAMEDICAL

SHINHWAMEDICAL

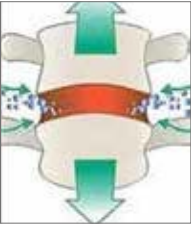


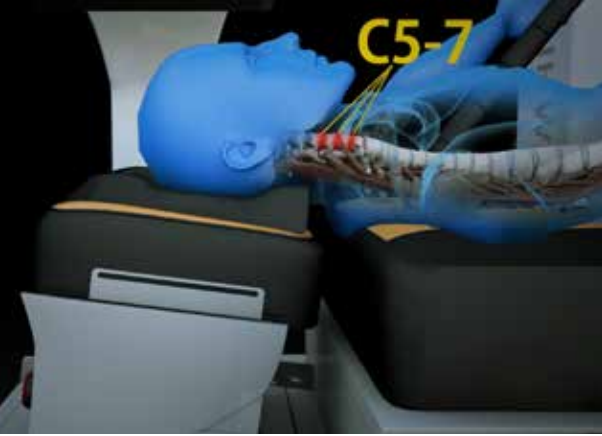
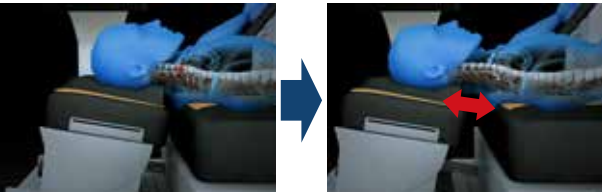
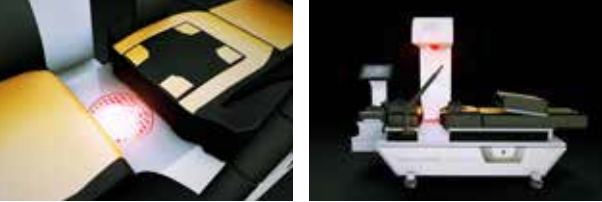
You can use it, depending on your purpose such as;

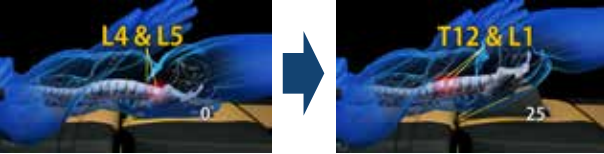
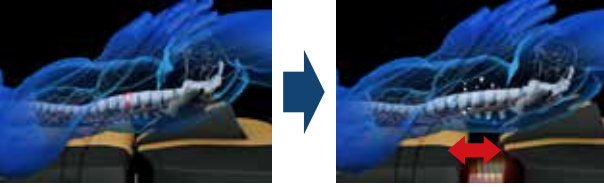
- 1. Basic chiropractic table + manipulation
(James M..Cox & James Cyriax Technics)
- 2. Premium auto decompression with programmed Technics
- 3. Combination of 1 and 2 in the above.



CORE THERAPY

	CORE 1 : DECOMPRESSION • Feedback System • Comfortable, Effective
	CORE 2 : CORRECTION • Target Set-up • Function : ①Unilateral ②Rotation ③Level
	CORE 3 : EXERCISE PRESCRIPTION • Intensity set-up : ①Low ②Medium ③High • Angle Set-up : 5 Steps
	CORE 4 : PROGRAM • Automatic Mode : Variously sophisticated Treatments
INDICATIONS • HNP & Degenerative Disc • Anti-colytic Spondylitis • Scoliosis & Stenosis • Joint Problems • Muscle Imbalance for rebuild	EFFECTS • Reduce pain • Increase vertebral flexibility • Reduce muscle spasm • Relieve nerve compression • Resorption of protruding discs

Part/Function	Image	Explanation
Cervical Target Angle		• Treatment principle (flexion-distraction) of Cox • Function to cure the cervical region to be treated through target angle adjustment • Adjustable from 0 ~ +15°
		• Functions according to medical examination and treatment principle • Functions to apply treatment according to the herniated direction of the disk • Up to 80mm adjustable for right & left
Cervical Decompression therapy		• Application of feedback system • Application of exercise function • Application of decompression function • Application of automatic program
Infrared Radiation therapy		• Up / Down Infrared ray is radiated to the rigid patient's body to relax the shortened muscles during the treatment. • The function maximizes the performance of the treatment. • Up infrared ray is optional.

Part/Function	Image	Explanation
Lumbar Target Angle		• Treatment principle (flexion-distraction) of Cox • Function to cure the lumbar region to be treated through target angle adjustment • Level adjustable from 0 ~ + 25°
Lumbar Unilateral application		• Functions according to medical examination and treatment principle • Functions to apply treatment according to the herniated direction of the disc • Up to 90mm adjustable for right & left
Lumbar Rotation application		• Treatment principle of Siriaks (orthopedic surgeon) • Treatment according to the type of escape of the disk • Adjust the angle of lumbar motion • Stop position from 0 to 18 ° for twisted spine & pelvis
Lumbar Decompression therapy		• Application of feedback system • Application of exercise function • Application of decompression function • Application of automatic program