


CHIROPRACTIC TABLE - 3 DROP MANUAL

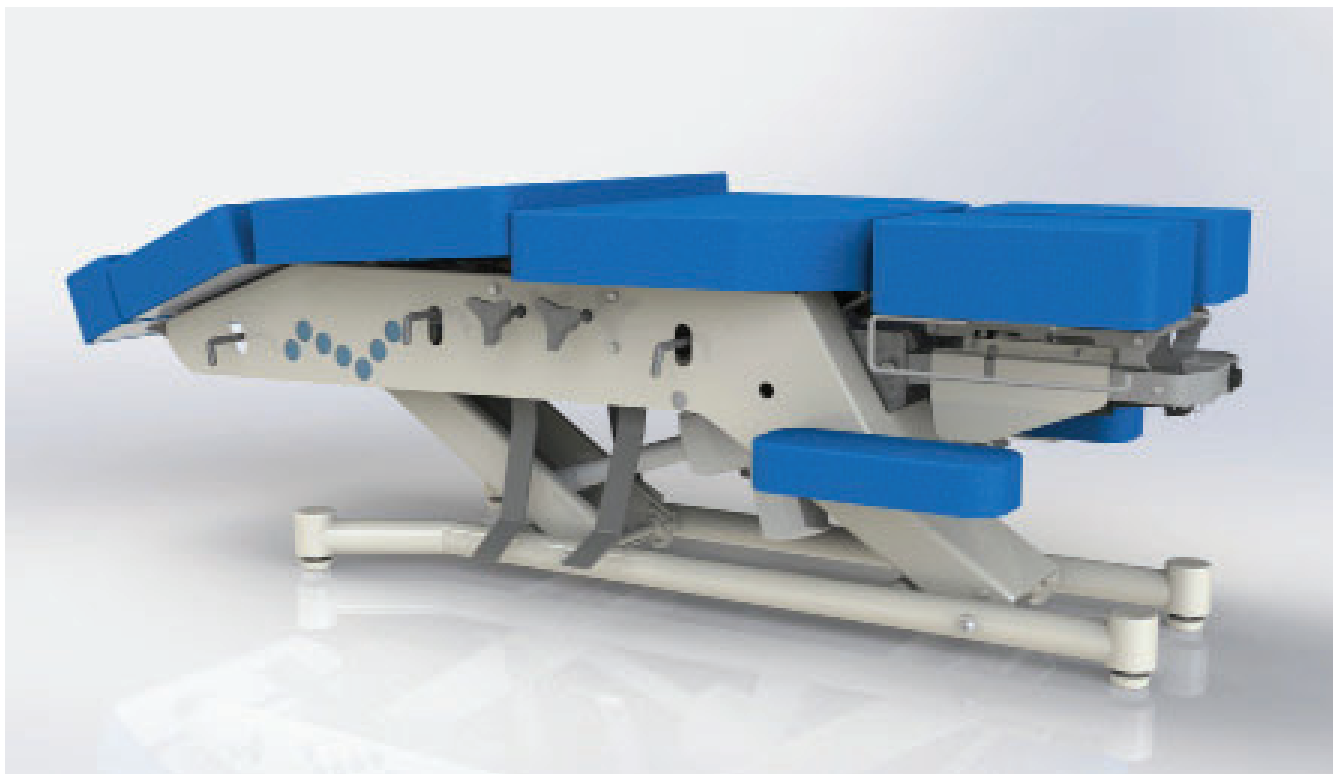

BODY DROPS

GV - II can be configured with up to three drops distributed by the main body. These are thoracic, lumbar and pelvic. In the stretchers with three and four drops there are two types depending on their movement on the one hand the thoracic and lumbar drops produce a vertical movement while the pelvic allows a regulation of the positive inclination angle. The actuation of these drops is carried out in a very comfortable way by the foot on the part of the user. The pedals for loading each drop are located on both sides of the stretcher for greater user comfort. The activation of the drops can be done in any position of inclination of the pelvic or dorsal cushion when you have that option. All drops have tension adjustment on both sides of the stretcher for greater comfort. The number of drops to be included in the stretcher will be optional by the client when placing the order.

MALFUNCTION HEAD

The head has a manually operates toggle drop with hardness adjustment using a knob on the left side. In any of the angles in which the head is positioned, the drop can be activated. For patient comfort, the separate on of the head support cushions is adjustable to adapt it to the patient's face. The multifunction head has a position for the paper roll and a guillotine at the opposite end for a homogeneous cut. Guidance of the paper performed with a central rod between the two cushions.

ORDERING DETAILS

CODE:	DESCRIPTION
GVII 3D	Chiro Practic Manual Table

