

VOGUE 7.0 TECAR

Fisiowarm 7.0 is more than just an extraordinary device. It's everything that you have ever wanted in a tecar therapy unit and didn't think it was possible. All this with a beautiful design that makes the difference. Fisiowarm 7.0 ensures the ultimate tecar therapy experience because we believe that it can't get much than this.

THERAPEUTIC INDICATIONS

Heat provides relief in chronic traumatic and inflammatory conditions, increases blood flow by encouraging cellular nourishment and the disposal of toxins, and improves the elasticity of connective tissue. The vasodilation obtained with the use of heat leads to an increase in the exchange of substances that promote the drainage of the inflamed site and the local influx of cells responsible for restorative processes.


- Contusive or distortive trauma
- Post-surgical rehabilitation
- Sports traumatology
- Muscle and tendon pathologies
- Scarring outcomes
- Pathologies of the spine and peripheral joints
- Pathologies of the vascular and lymphatic system
- Pathologies of the nervous system of orthopaedic origin
- Rehabilitation of the pelvic floor



EFFECTIVE GEOMETRIES OF TECAR THERAPY

By Geometry in Tecar therapy, we mean the path the current takes in the patient's body with respect to the position between the plate and the electrode which are those that see the parallel positioning between plate and electrode, where the geometric projection of the electrode ends in the centre of the return plate.

With this type of geometry the passage of the current within the body and how the energy will be yielded is very clear. In fact, in a treatment using the resistive system in correct geometry, we will only have to know the anatomy of the structures that are between the plate and electrode, and according to the different features of the tissues we will know where and how the energy will be yielded.

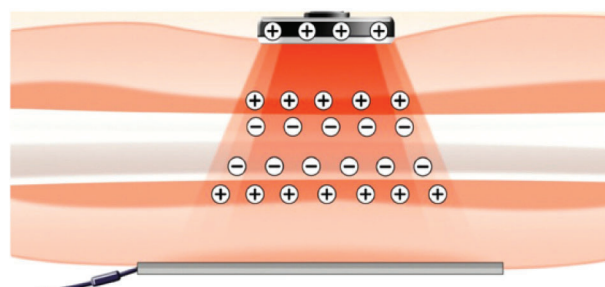

Schema di trattamento in
'Corretta Geometria'


ORDERING DETAILS

CODE:	DESCRIPTION
VOGUE 7	Tecar Therapy With Ipad

UNIQ: CAPACITIVE + RESISTIVE

Uniq combines both the capacitive and resistive treatment methods in a single unipolar electrode.

It enables same-session treatment both of surface tissues (capacitive electrode), and deep tissues (resistive electrode), halving treatment times and at the same time doubling efficiency.

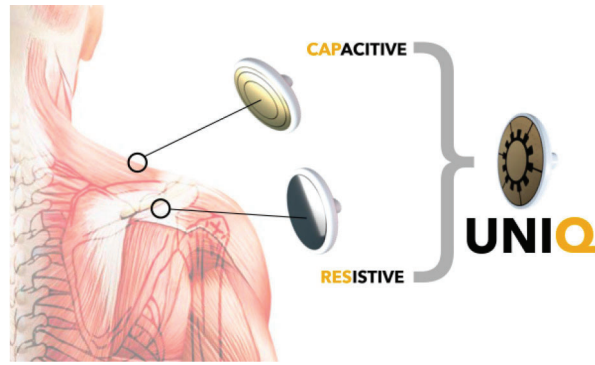

Tecar Therapy WITH 28 Accessories to cover the full range of applications (musculoskeletal, Aesthetic and Pelvic floor disorders)

- Unipolar
- Bipolar
- Myo facial
- Intimity
- Physio Aesthetic
- Neurodynamic
- Trigger point

VERSION FEATURE



Power	300W	400W
OPERATOR INTERFACE	Soft-touch resistive panel 10.1 inch tablet	Soft-touch resistive panel 10.1 inch tablet
OPERATING FREQUENCY	0.5/0.75/1MHz	0.5/0.6/0.7/0.75/0.8/0.9/1MHz
OUT PUT	Continuous 50% continuous Pulsed 1-1	Continuous 50% continuous Pulsed 1-1 Pulsed 1-2 Pulsed 2-1
OPERATING MODE	Unipolar Bipolar Neuro Dinamico Fixed Automatic	Unipolar Bipolar Neuro Dinamico Fixed Automatic
ACCESSORIES PROVIDED	1 hand piece for res/cap unipolar apps 3 resistive unipolar electrodes 3 capacitive unipolar electrodes 1 hand piece for res/cap bipolar apps 3 capacitive bipolar electrodes 1 resistive bipolar fixed electrode 1 neutral plate 1 power cable 1 Plexiglas accessories holder	1 hand piece for res/cap unipolar apps 3 resistive unipolar electrodes 3 capacitive unipolar electrodes 1 hand piece for res/cap bipolar apps 3 capacitive bipolar electrodes 1 resistive bipolar fixed electrode 1 neutral plate 1 power cable 1 Trolley with built-in accessories holder.
COMPATIBLE ACCESSORIES	All	All
TRANSFER	Capacitive Resistive	Capacitive Resistive
PRE SET PROGRAMME	22 Operating Protocols	46 Operating Protocols
SERVER SIDE SERVICES	Included	Included
DIMENSIONS	W32 x D30 x H11	W32 x D30 x H11
TROLLEY DIMENSIONS	W56 x D x 56 x H140	W56 x D x 56 x H140
WEIGHT	4.5 kg	4.5 kg